

User Manual



HEALTH • TECHNOLOGY • ECO-FRIENDLY



Contents

Safety Instruction	1-2
User Guide	2-3
Warranty	3
Features and Usage Demonstration	4
Operation and Instructions	5-6
Names of Parts	7
Product Maintenance	8
Troubleshooting	9
Product Specifications	10

Thank you for purchasing this Leg & Foot Massager. Please read the operation manual thoroughly for the correct methods of usage and take special notice of the safety section before you start. This manual should be kept on hand for easy reference on safe and efficient operation of the massager.

■ Safety Instruction

Important Safety Instruction

1. To ensure safe and correct use of the massager, do not operate it without carefully reading this operation manual.
2. Close supervision is necessary when the product is used by any person who is an invalid, handicapped, or with any form of physical, sensory or mental incapability and with limited experience and knowledge of the product(including children).
3. This product is not intended for children, ensure that they do not play with the product.
4. There are heat generation components inside, the surface of the product could be hot, so the people who is not sensitive to heat should be more careful when using this product.
5. If the wire is damaged, to avoid danger, one should ask manufacturer or authorized maintenance institution/team to change the wire.
6. The massager should be kept away from sharp items. Do not use the massager when the fabric cover is torn in order to prevent toes be drawn into the machine causing injuries. One should ask manufacturer or authorized maintenance institution/team to change or repair before use again.
7. Please unplug the machine before cleaning or after finish using.
8. Do not cover the massager during operation. This may result in a rise in motor temperature and subsequent malfunction or damage.
9. Stop using immediately should you experience any discomfort.
10. When the machine malfunctions, please send it back to manufacturer or authorized maintenance department to repair before you use it again.
11. Please do not use the spare parts unless it's recommended by original manufacturer.
12. Never drop or insert any object into any opening of the product.
13. Do not operate the massager with wet hands.
14. Do not expose the control switch to water or any liquids.

15. To disconnect, switch off all controls and remove plug from electrical outlet.
16. This product is not intended for self-diagnosis or self-treatment for conditions that should be managed by a qualified health care provider.

Usage environment

1. Do not use the massager in an excessively dusty or humid environment such as a bathroom since this may result in malfunction, electric shock or damage to the product.
2. Do not use the massager in a room with temperature of 40°C or higher. Do not expose the massager to heaters, stoves or direct sunlight.
3. Do not use the product when the ambient temperature is below 5°C, for the thickening of the lubricant could increase abrasive wear of the gears, which will shorten the operational life span of the product.
4. Better not use the product outdoor or place it under direct sunlight, it might cause damage to the paintcoat of the product.
5. When moving the machine from low temp environment to relative high temp to use, the inside metal parts might malfunction or even damaged because the water vaporization. We suggest placing the machine under room temperature for more than one hour before power it on.

Following person please consult with doctors before use

1. Patients that are under treatment or post-surgery patients.
2. Patients who suffer from high blood pressure, cardiac disease, cerebral arteriosclerosis or cerebral thrombosis.
3. Patients who suffer from malignant tumor
4. Pregnant women or women during menstruation.
5. Patients who suffer from osteoporosis.
6. Person with cardiac pacemakers or have implantable medical device.
7. Patients who suffer from acute diseases.

■ User Guide

Use this massager 1-2 times a day, 5-10 minutes each time would effectively eliminate muscle fatigue, improve blood circulation. Use this massager for 5-15minutes before sleep, would effectively improve sleep quality.

Notice

1. Do not fall asleep during using.
2. Please do not use this product when you are inebriated or when you are not feeling well.
3. Please do not massage over 30mins, warm foot do not over 2 hours.
4. If the cloth cover is broken, please don't use it.
5. Do not use this product one hour before or after meal, it might cause opposite effect.
6. When moving the machine from low temp environment to relative high temp to use, the inside metal parts might malfunction or even damaged because the water vaporization. We suggest placing the machine under room temperature for more than one hour before power it on.

■ Warranty


1. The right of warranty is not transference.
2. The warranty of this product is one year, starting from the date user purchase the product. During warranty, maintenance is free based on valid receipts(the appearance is not within warranty, if need mailing, users only need to pay for the shipping fee).
3. The warranty of this product is under normal use for a period of ONE YEAR.
4. Original purchase receipt and warranty card is required when claiming warranty.
5. The purchase time is subject to the valid documents issued by authorized seller.
6. This product is only suitable for household use, renting or use for business purpose are not within warranty.
7. The following conditions do not belong to the scope of the warranty, only to provide paid services (charge labor costs and material costs)
 - a. All damages caused by human factors, including being used under improper working condition, used and kept not according to the manual, and damages caused by improper use etc.
 - b. Disassemble, repair, modify the massager by users or by non- authorized third party.
 - c. Damage caused by bad transportation.
 - d. Damage caused by other irresistible factors (such as floods, fires, lightning and earthquakes etc.).

■ Features and Usage Demonstration

Features

1. This product is developed by professional product R&D team. The design is elegant and with assorted colors. Classy lacquer surface radiates fashionable qualities. It is also light-weighted and portable.
2. Ergonomic design, enable you to enjoy professional leg & foot massage at any time. The ergonomic design, unique styling and acupoint distribution all fit to the human leg & foot structure. Just like bringing a professional massagist home, you will enjoy ultimate relaxed massage on your leg & foot while watching TV program, or before sleep.
3. Double safety protection makes you feel secure when use this product. When the product is overheating or the current is overload, the circuit will break off automatically. This double protection is safe and reliable, ensure you to enjoy the massage joyous and carefree.
4. Combined various massage ways of kneading, scrapping, shiatsu, airbag press, warm foot etc.
5. Unique calf massage and care function.
6. Nano removable washable foot cloth cover, keep you clean at all time.
7. Equipped with double sensitive touch panel.

Usage Demonstration

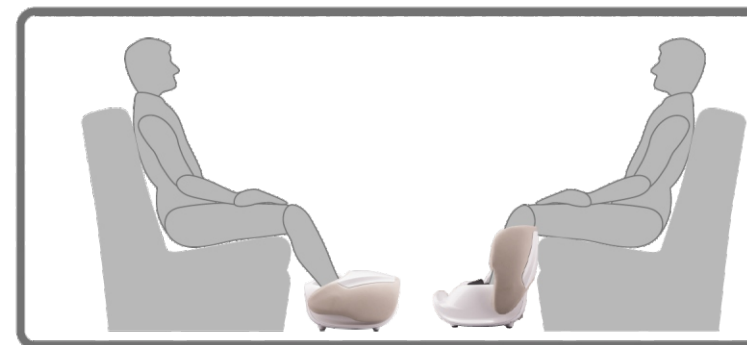

■ Operation and Instructions

Operation



ON/OFF :

Press the ON/OFF button on the front panel of the massager to start the massager and enter the automatic working state.



HEATING :

Press the HEATING button to turn off or turn on the foot heating function.



STRENGTH :

Press the STRENGTH button, you can choose the appropriate air massage strength, the button indicator light according to the weak, medium and strong, synchronized display for the blue, blue and red, red; long touch the button for one second, you can turn off the air massage function.



SCRAPING :

Press the SCRAPING button, you can turn off or turn on the foot massage function; long touch the button for one second, you can turn off the foot massage function and switch to a soft foot air care experience.



VIBRATION :

Press the VIBRATION button, you can turn off or turn on the foot vibration function and select the vibration education experience for different strength modes.



LEGS :

Press the LEGS button, you can turn off or turn on the leg foot massage function. Long touch the button for one second, you can turn off the foot massage function and switch to soft leg foot air care experience (This function is not available when the clamshell is folded) .

Instructions

1. Place the massager on a clean floor and check if the machine's foot cloth cover are loose or broken
2. First plug the end of the adapter into the massager outlet, the other end connected to the power cord, and then plug the power cord into the mains socket, the massager power indicator light starts flashing, the massager into the standby state, people sitting in the rear of the massager, put the feet into the product foot cloth cover (note not stand to use).
3. According to the needs of massage, you can choose to flip part fold or unfold, flip folding can massage the foot, flip unfolding can massage the foot and calf at the same time.
4. The upper operation panel is invalid when the flip is folding, the lower operation panel is invalid when the flip is unfolding.
5. If you want to change the flip state in the using, you need to press the ON/OFF button, until the product stops working and then lift the flip.
6. In the working state press the power button or the end of the timer time, the massager will stop working, the power indicator flashes, enter the standby mode.




The massager is finished using in the clamshell state(or stop in the middle), should not be flip folding immediately, the massager should be allowed to stand for one minute, until the gas is completely released and then folded.

Tips:

After use, a cup of warm water is better.

To prevent children accidentally start the massager caused an accident or affected parts, please unplug the power cord from the electrical outlet when you don't use the massager, and proper storing.

■ Names of Parts

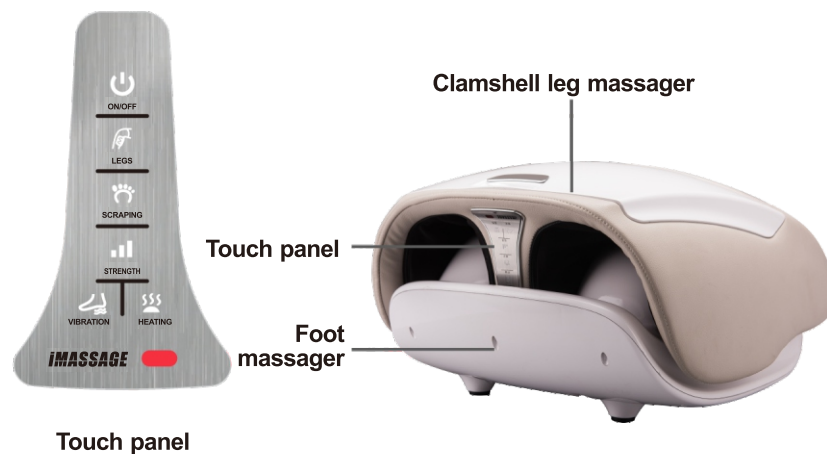

■ Product maintenance

Storage

1. When not in use, place the massager in packing box and place it in a safe, dry, ventilated, cool place.
2. Do not twist the wires.
3. Do not use wire hanging massager.
4. Avoid sharp objects scratching the massager's surface
5. Do not place the massager in direct sunlight or high temperature.

Clean

1. Before cleaning, make sure the power switch is off and unplug the power cord.
2. Please use a dry cloth to clean the surface of the massager. Do not use corrosive cleaners, such as benzene or thinner to clean the massager.

Maintenance

1. Keep the massager clean and away from high temperature and humidity place.
2. If there is a long time not to use the massager, please use the dust cover (please buy their own) trap to prevent dust from falling into.



This massager has a temperature protection device, it will automatically stop working when the temperature exceeds the safe temperature, please unplug the power to allow its natural cooling and then use.

■ Troubleshooting

If you have problems using this massager, please refer to the following tips. If this massager is still in a failed state, return the massager to the designated service unit or an authorized technician for inspection and repair

Abnormal Situation	Possible Causes	Treatment Methods
Not working	1. Not power on 2. The function button is not turned on	1. Connect the power supply or check whether the power supply is connected properly 2. Turn on the ON/OFF button
Suddenly stop operating	1. Power off 2. Massage time more than 15 minutes 3. Continuous using time is too long, massager overheating protection 4. Force too much, massage ball stalled	1. Check if the power connector is connected 2. Turn on the ON/OFF button again 3. Turn on the On / Off button again after cooling the massager 4. Relax the feet, reduce the load on the massager

■ Product Specifications

Product Name	Leg Massager
Rated Voltage	DC24V
Rated Power	55W
Timing time	15mins

- ▲ Our company will continue to improve product performance and appearance, subject to change without notice.

Special statement

1. All contents of this manual have been carefully checked, if any print errors or misunderstandings on the content, our company reserves the right to explain
2. If the product has technical improvement, will be compiled into the new version of the manual without notice, product appearance or color changes, please prevail in kind.